

# TECNO-GAZ

## NEWS

Number 2

Information and opportunities  
directly from the  
manufacturer to the user

EXTEND YOUR VISION. IMPROVE YOUR WORK WITH

# TECNO-GAZ PROJECTS

TECHNOLOGY  
BUILT  
AROUND  
YOU.

SURGICAL  
CLINIC  
PROJECT

BIOLOGICAL  
RISK  
PROJECT

CONSCIOUS  
SEDATION  
PROJECT

FIRST AID  
PROJECT

FOCUSON

### The New MultiSteril Fast

Faster, safer and cleaner than ever before. Discover how to make the most of your time with this Tecno-Gaz unique device.



Page 3

### KYRI surgical Chair

A revolutionary approach to dental surgery



Page 5

### Video: Sedative analgesia

#### Learn from experts

Check our new playlist with a complete overview of everything you have to know



Page 12



TECNO-GAZ  
INTEGRATED PROJECT

# A WORLD BUILT AROUND YOU.

We offer a *complete Tecno-Gaz clinic.*

Tecno-Gaz offers a wide a range of solutions to handle, manage and expand the operations of your practice. Our solutions not only provide the single product, but also a range of services, skills and expertise.

# 360



# THE NEW MULTISTERIL FAST



With **Multisteril Fast**,  
sterilization in only **10'**



- + Exactly the time required to order and wrap your instrument sets
- + The presence of the assistant is only required for wrapping.
- + Multisteril Fast, totally automatic and autonomous, will carry out all the pre-disinfection phases by itself.
- + Ask your dealer for a presentation, thus saving precious time for your patients.

**Faster** ✓  
**Safer** ✓  
**Cleaner** ✓



TECNO-GAZ  
STERILISATION  
PROJECT

# Sterilisation FAST



Complete  
Cycle

25'

NEW!

**MULTISTERIL FAST**  
COMPLETE CYCLE  
in 25 minutes!

1

**MultiSteril  
Fast**



**MULTISTERIL CD**

Liquid for disinfection and cleaning



25'

COMPLETE CYCLE

STAND-BY  
DISINFECTION

NEW!

ULTRASOUND  
CLEANING

AUTOMATIC RINSING

DRYING

**Multisteril + Multisteril CD**


All phases are automatic: disinfection, cleaning, rinsing and drying. Safety, time saving, money saving.

The first Certified System:

**Certified cleaning**

**Certified disinfection**

**Instrument protection**

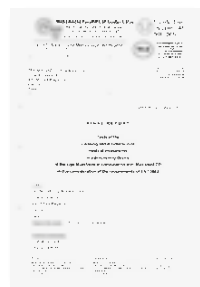

# ADOPT TECNO-GAZ PROTOCOL!



WITH ONLY 3 PIECES  
OF EQUIPMENT  
YOU CAN HAVE THE  
BEST STERILISATION  
PROCESS!

MAXIMUM  
OPTIMIZATION  
AND STREAMLINING  
OF THE LOADS

1 FULL CHARGE OF  
MULTISTERIL UNIT

IS THE SAME AS 1  
FULL CHARGE OF  
AUTOCLAVE



(1) **MULTIFUNCTIONAL UNIT  
MULTISTERIL**

(2) **AUTOMATIC THERMOSEALING  
MACHINE  
TECNOSEAL**

(3) **AUTOCLAVE B TYPE LOADS  
ONYXB**

Time saving

**+50%**



**TecnoSeal**

**Onyx B**

Available in 5 and 7 kgs

Patented P.I.D. automatic  
adjustment system

SD Card for data tracing and  
autoclave upgrades

User validation at cycle end

Baric alignment

Autoclave with large capacity

Fast cycles and small loads cycles

Best available components

Simple to use

"Full Load"  
System

7 Automatic  
Cycles

SD card

PID System



Automatic drive system  
Photocell (10 m./Min.)

Maximum productivity,  
minimum time

Electronic board  
Management control for  
use parameters





# KYRI

The Kyri project can be configured according to a purely surgical logic, which can then be developed in a multifunctional surgery, with the addition of the water unit and the cart with instruments.



## KYRI

Kyri is the only surgical chair designed and built according to the rules of dental ergonomy. The chair movements are controlled by four motors, thus allowing the chair to be moved into any position, facilitating the operator's work and increasing the patient's comfort. The chair can be equipped with a satellite arm for cannula suction and side rails on which it is possible to put any accessory: armrests, infusion supports, supports for housing arms, etc..

### 3-axis orbital headrest

70% of implants involving molars and premolars require lateral surgery. This headrest ensures better posture and less muscular pressure for the patient, thereby facilitating such surgery.



### Infinite number of positions based on 5 axes

Kyri assists you in any position. Thanks to its 5-axis movement, it can assume any position required to ensure its conventional operation, as well as unusual or unexpected positions.



Dental

### Anti-shock and emergency positions

Orthopneic or semi-seated position

Trendelenburg

Fowler

### Surgical Table

Variable height

Multi-purpose  
Surgical Clinic

Surgical Clinic

# A modularity that allows you to expand your surgery according to your requirements.

## + Ergoseat

The device is designed to guarantee correct posture, thanks to the inductive system with integrated motion between seat and back.



## + Smart Trolley

The first range of technical furniture, providing greater efficiency of your daily work. A real innovation and evolution in dental cabinets.



80+

More than  
80 possible  
configurations



## + Orion 40Ds

Universal lamp, with field selection: 15 and 20 cm field width adjustable light intensity, intense spot light depth, double colour temperature, extreme manageability.



## + KYRI CART

Possibility of configuration and operational flexibility thanks to Kyri Cart. Five slots for instruments and total control of Kyri settings to store and recall your preferred positions.



## + Kyri Water Unit

The element that makes the surgery truly versatile, (from surgery to general dentistry). It can be placed under the worktops and equipped with suction cannulas.



## + Kyri Dss

New surgical aspirator: adjustable suction, control pedal, double basin and disposable bags with valves.



## + Kyri Bridge

The bridge table is a simple accessory positioned in the common ergonomic space between doctor and assistant. Its strategic position contributes to the quality of surgical treatment.



**80+**  
More than  
80 possible  
configurations

# SMART TROLLEY

Unlimited  
setups



Tecno-Gaz has developed the SMART TROLLEY concept, changing what seemed unchangeable. SMART TROLLEY is a real revolution, not only because it completely innovates a sector, but also changes the identity and role of a major component of the dental surgery. The 'FURNITURE' is no longer simply a an attractive element, but a component of organisation, order and efficiency - an asset in the daily work of the surgery.

**IT IS TIME  
TO MAKE  
YOUR CLINIC  
OPERATIONAL**

- ✓ **Created to make work easy and productive**
- ✓ **From the container to the operating device**
- ✓ **Suitable for all clinics: can be easily integrated and configured**

The first impressions and moments of any treatment are the most important factors for both patient and operator. This is when levels are set for tension and the relative quality of any following dental procedure. A clinically welcoming environment provides the patient with a positive impression and the dentist with excellent ergonomic operation.




## + SMART TROLLEY Intra Operator



Plan in case Surgery does not go as planned. You may need different implants size, additional instruments, anesthetics, ... Because of surgical time to respect, hygiene procedure, and patient management, it is not thinkable to leave the room to look for the missing components. Smart Trolley Intra Operator is the temporary stock you need, at immediate reach.

HQ guides with  
removable  
drawer

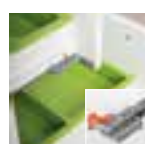

## + SMART TROLLEY Electromedical



Tilting  
surfaces

Dental staff use a large amount of electromedical instruments during dental procedures. The solution is the Electromedical Smart Trolley which, with its centralised power supply system, modularity, and the accommodation of bulky pedal commands, becomes an indispensable organisational and functional element.

## + SMART TROLLEY Postoperative



To save time and efficiently manage waste and contaminated material, we have developed the Postoperative Smart Trolley. Thanks to its built-in decontamination tank and recycling compartments, it is a hub for the collection and sorting of infected material.





Light intensity up to:  
**130,000 Lux**

Colour temperature:  
**4,500/5,000 K**

Illuminated field:  
**15 cm and 20 cm**

Special for dental use

## Why use a dental surgical light?


The lighting must provide precise transmission of visual information from the work area to the surgeon during the operation, but must also guarantee eye protection.

## Possible visual damage by traditional lighting

Light is essential for the eyesight, which must be taken care of from the beginning of your career, not when problems arise.



Invest in your main diagnostic tool:  
**your eyes.**

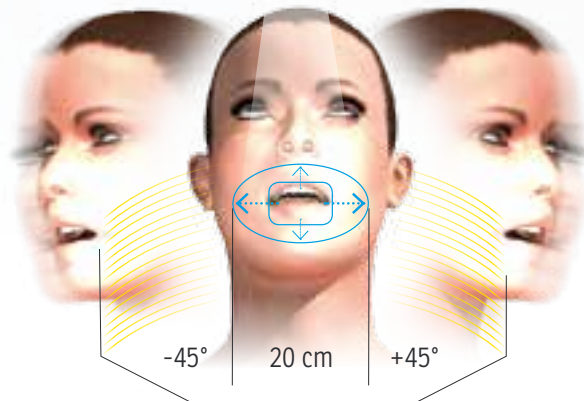

# PROTECT YOUR EYES

with a dental surgical light



## The unique lamp specifically designed for dentistry!

Orion 40 DS is a development of the famous Orion 40 dedicated to the clinical needs of the dentist. The first surgical light specifically developed for dental use thanks to photometric measurements calibrated according to the dentist's requirements: 80 cm focus, 15 cm illuminated field up to 130,000 lux and temperature of 4,500/5,000 k. ORION 40 DS is the first multifunctional light for dental practices, which, thanks to the unique E-Dentalight technology (joint function), switches to surgical light.

## E-Dentalight

This function enables the replacement of lighting generated by traditional lamps. This unique feature enables the activation of two modules and the adjustment of their intensity and colour temperature.

## Cold light

The increase of heat radiation on the surgical field is proportional to the quantity and power of the LEDs used. Thanks to the indirect light concept, Tecno-Gaz uses fewer LEDs.

## Parabolic indirect lighting

## Double illuminated range

The ranges of light intensity of the Orion 40DS Lamp change according to operator's requirements:

- E-Dental Light Function: 15 cm field with variable intensity from 20.000 to 55.000 lux for general dentistry treatments.
- ORION Light Function: 20 cm extended field with variable intensity from 40.000 to 130.000 lux for surgical treatments and implants placing.

## ORION40 DS

Free standing  
Art. LC103LRD

Wall  
Art. LC101LRD

Ceiling  
Art. LC102LRD

Ceiling mount for camera  
Art. LC106LRD

## ORION60

Light intensity up to:  
**160,000 Lux**

Colour temperature:  
**4,500/5,000 K**



Reflector external diameter:  
**Ø 630 mm**

## ORION80


Light intensity up to:  
**160,000 Lux**

Colour temperature:  
**4,500/5,000 K**

Reflector external diameter:  
**Ø 790 mm**



**SEDATIVE  
ANALGESIA**

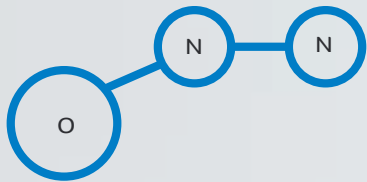

More value with

**Service Plus**

**YouTube Video available**

Page 12

**NEW**

# MasterFlux Plus

## New IntelliFlux delivery system

IntelliFlux is designed to solve all the problems caused by sedation mask circuits. IntelliFlux has a connection with backflow valve to prevent dispersion and increase and enhance the sedation process.



**IntelliFlux**  
DELIVERY SYSTEM

NEW INTEGRATED SYSTEM  
FOR THE SUPPLY OF GAS

ITS SPECIAL DESIGN PREVENTS A NEGATIVE  
VISUAL AND PSYCHOLOGICAL IMPACT ON THE  
PATIENT

NEW CIRCUIT: BETTER COMFORT  
FOR PATIENT AND OPERATOR



## Automatic flow control

The system has been designed to make flow administration quick and easy. Turning the O<sub>2</sub> knob, the flow required by the patient is set in litres per minute, depending on his/her respiratory volume. The special regulation knob enables the percentage of nitrous oxide delivery to be automatically adjusted from 0 to 70%.

The columns indicating the flows have been designed for fast and easy reading. The couplings have been fitted on the back of the flow control box and their attachments are physically different to safeguard against the risk of inverting the hoses for the two gases. Yet another precaution to guarantee safety in the use of the machine.

## Safety device

## Minimum oxygen concentration

A constant percentage of 30% oxygen safeguards against the risk of HYPOXIA and ANOXIA.

## O<sub>2</sub> Emergency

By pressing the O<sub>2</sub> device, thus cutting of the N<sub>2</sub>O flow, **100% oxygen is delivered both to the control balloon and to the patient being treated.**

## Anti-expired gases valve

A special valve prevents expired gases from being recycled in the circuit. Additionally, a special device enables the patient to breathe ambient air in the event of a gas failure.

## N<sub>2</sub>O cut-off

Cuts off the N<sub>2</sub>O flow when O<sub>2</sub> pressure is zero (no oxygen).



# THE NEW CIRCUIT

IntelliFlux supplies gas through two central tubes, which do not obstruct operations.

Rear elastic support to secure the mask to the patient's face.

## Ergonomics

- Perfectly suited for 4-hand usage.
- Circuits joint located above the patient's head avoiding impediments to the treatment.
- Can be moved in all positions.

## Fast and efficient sedation

- Perfectly sealing mask suitable for any face shapes and curves for a faster, more efficient sedation treatment.

## Pleasant treatment

- The patient is sedated, enjoying a peaceful experience, breathing through a mask that will leave no marks on his/her face.
- The frontal holder and the rubber elastic band ensure stability.

## No gas dispersion in dental clinic environment

Many particular specifications contribute to prevent dispersion such as:

- Backflow valve.
- Adaptable connector to dental suction unit.
- New easy-breathe membrane.
- Face-adaptive, single use silicone form.

## Advantages of the new mask and new circuits

Nose masks are no longer 100% silicone, but have a rigid transparent polysulfone centre shell, which enables visibility and prevents any torsion, thus guaranteeing a perfect fit. The adaptive cushioning and front forks enable perfect, stable positioning, thus preventing any gas dispersion or leakage. All parts of the nose mask are autoclavable and the centre connection contains the backflow valves.

STANDARD



1

STANDARD



2

STANDARD



3

OPTIONAL



4

- Gas dispersion control in the clinic environment
- Perfect sealing on any face
- Greater comfort for patient and operator



NEW!

Pediatric MASKS

NEW!

size 1



NEW!

size 2



# THE NEW MASKS

**SEDATIVE  
ANALGAESIA**

# Services and Resources

A WIDE RANGE  
OF **EXCLUSIVE**  
SUPPORT SERVICES  
TO **IMPROVE YOUR  
WORK**



## Free Video Tutorials online

Tecno-Gaz has made 14 videos, describing the main advantages of the daily use of sedation. These are available on the TECNÓ-GAZ YouTube channel. Also on demand by writing to [marketing@tecnogaz.com](mailto:marketing@tecnogaz.com)

- Individualised dosage - Sedation level & incremental technique - Induction and recovery - Safety margins
- Metabolism - Side effects and contraindications - Analgesic potential - Psychically and physically disabled patients - Children - Machine safety - Treatment safety - Long, complex, traumatic treatments - Fear, phobia, anxiety

## Service Plus for Sedative Analgaesia

- 1 In the countries in which Master Flux is distributed, we can count on the collaboration on the part of universities and professionals, able to offer information and training on the use of our products, to implement the technique in daily application.
- 2 In Spain, the collaboration with UIC (Universidad Internacional de Catalunya) is very fruitful and, every year, this university proposes conscious sedation training courses.
- 3 The sedation system is an integrated part in the work of the university.
- 4 The professors of the medical faculty of the university collaborate actively with our company by holding conferences and presentations regarding conscious sedation: Manuel Ribera, Victor Gil, Inma Puga etc.



## Conferences

Every year, Tecno-Gaz organizes conferences about the advantages of the use of sedation in dentistry, in partnership with our distributors.



## Tutoring

Tecno-Gaz and its distributors operate on the territory with qualified representatives, who can provide continuous on-site training on conscious sedation and products. The distributors can also have equipment in their showroom, where they can show the whole TECNÓ-GAZ product line.



## Technical support

Tecno-Gaz has a wide network of dealers, who can provide on-site technical assistance.



## End-user training

To optimise the use of sedation equipment at a daily level, it is important for a user to express doubts and enquiries to an experienced professional. Through our distributors, we can organize specific courses for small groups of users in the presence of a university professor providing practical training.



## Bibliography

An extensive scientific bibliography confirms the safety of the technique in dentistry, and we can make all this bibliographical documentation available.



## University training

Through its distributors, Tecno-Gaz can offer university training and award an attendance certificate.



## Remote phone assistance

To optimise the use of sedation equipment at a daily level, it is important for a user to express doubts and enquiries to an experienced professional. It is possible to book a date and time for a 1 hour consultation about sedation and technique.



### Conscious sedation in dentistry



Doctor Victor Gil the collaborator presents: Metabolism

TECNO-GAZ

**A GREAT  
EXPERIENCE  
BECOMES A  
GREAT VALUE FOR  
CUSTOMERS**



The life-saving device  
**FIRST-AID**

IS IT AN **OBLIGATION** OR NOT?

# Defibrillator

THE UNIQUE DEFIBRILLATOR WITH **EXCLUSIVE FEATURES!**

Obligations are often seen as "constraints", discussed and interpreted according to one's own interests. But when we speak about life-saving products, we must focus on ethics and common sense, as opposed to the concept of obligation. It is clear that in all environments where surgical techniques are practised, as in the dental field, devices must be adaptable to intervene in case of critical involving the patient. Therefore, the pulse oximeter, oxygen and AED defibrillator, become instruments which must be present.

✓  
**UNIVERSAL ELECTRODES**  
adult/child

- 1 **Cover:** the cover is designed to protect the action icons, patient mode button and charge button.
- 2 **Status indicator:** the status indicator indicates the status of the unit, the temperature and the battery level.
- 3 **Paddle connector:** the paddle connector is used to connect the electrodes.
- 4 **Action Icons:** the LED indicator flashes red under the corresponding action icon.
- 5 **Patient mode switch:** once the user has identified the patient based on type, patient mode (either adult or pediatric) is selected using the mode switch.
- 6 **Charge button:** once preparation for the electric discharge is completed, the current discharge button will flash. Press the current discharge button to deliver the shock.
- 7 **Sliding button:** the sliding button opens the cover, to turn on push the sliding button to the right.
- 8 **SD card slot:** the SD card is required to save data and update AED firmware.
- 9 **Infrared communication port:** the infrared communication port is used to communicate with the PC.

## + Tecnoheart Plus

### SEMI-AUTOMATIC DEFIBRILLATOR

Tecnoheart Plus enables safe defibrillation of the heart muscle walls, as it has sensors which can recognise cardiac arrest due to arrhythmia, ventricular fibrillation and ventricular tachycardia.




## + E Point

Mobile First Aid Centre



## Vital Test +

Professional pulse oximeter designed for use in dental clinics.

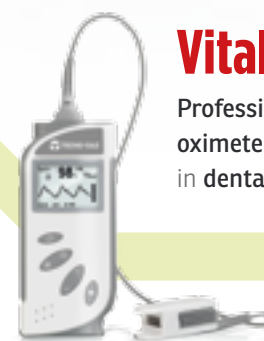

## + Air-Ox

Resuscitator



Web: [www.tecnogaz.com](http://www.tecnogaz.com)





# Welcome to the *Savings Area!*

## BUNDLE 1

**Multisteril Fast+  
TecnoSeal+  
Onyx B 5.0+Clean**

List Price ~~9.443€~~

Offer Price

**7.500€**



Art. AB015ZXY

**Onyx B 5.0**

Autoclave for the sterilisation of type B loads



Art. VM013ZVM

**Multisteril Fast**

Washer disinfectant



Art. TA007TAU

**TecnoSeal**

Automatic sealing machine



Art. 2070S

**Clean**

Smart manual instrument cleaning drum

## BUNDLE 2

**Optical data reader+  
Tecno-Print+Pura Plus+  
TecnoSeal+Onyx B 7.0**

List Price ~~8.332€~~

Offer Price

**6.800€**



Art. AB012ZXZ

**Onyx B 7.0**

Autoclave for the sterilisation of type B loads



Art. OI002NPU

**Pura Plus**

Reverse osmosis system



Art. TA007TAU

**TecnoSeal**

Automatic sealing machine



Art. 2110S

**Tecno-Print**

Printer for autoclave



Art. 3ELNL0001

**Optical data reader**





Art. VM013ZVM  
**Multisteril Fast**  
Washer disinfectant

Art. TA007TAU  
**TecnoSeal**  
Automatic sealing machine

Art. 2070S  
**Clean**  
Smart manual instrument cleaning drum

**BUNDLE 3**

**Multisteril Fast + TecnoSeal + Clean =**

List Price ~~5.466€~~

Offer Price

**3.900€**



Art. AE005Z12  
**TecnoHeart Plus**  
Semiautomatic defibrillator

Art. 1600AS  
**MasterFlux Plus**  
Sedation system with manual fluxmeter box

Art. 11250  
**Vital Test**  
Pulse oximeter

**BUNDLE 4**

**MasterFlux Plus + Vital Test + TecnoHeart Plus =**

List Price ~~6.293€~~

Offer Price

**5.500€**



Art. MF835ZMF / MF836ZMF / MF837ZMF / MF820ZMF / MF824ZMF

**Extra set of Masks and Circuit**

Art. 1600AS

**MasterFlux Plus**

Sedation system with manual fluxmeter box

Art. 11250

**Vital Test**

Pulse oximeter

**BUNDLE 5**

**MasterFlux Plus + Vital Test + Extra set of Masks and Circuit**

List Price ~~5.457€~~

Offer Price

**4.500€**



Art. MF835ZMF / MF836ZMF / MF837ZMF / MF820ZMF / MF824ZMF

**Extra set of Masks and Circuit**

Art. MF001EMF

**Master Flux Electronic**

Sedation system with digital control fluxmeter box

Art. 11250

**Vital Test**

Pulse oximeter

**BUNDLE 6**

**MasterFlux Plus Electronic + Vital Test + Extra set of Masks and Circuit**

List Price ~~7.257€~~

Offer Price

**5.500€**





Art. M0002MPS  
**E-Point**  
Emergency cart

Art. AE005Z12  
**TecnoHeart**  
Automatic defibrillator

## BUNDLE 7

**E-Point + Air-Ox +  
TecnoHeart Plus +  
Monitor MD**

List Price ~~6.411€~~

Offer Price

**4.500€**



Art. 3002-E  
**Air-Ox**  
Portable oxygen



Art. 930325  
**Monitor MD**  
Electronic  
multiparameter monitor



Art. AE005Z12  
**TecnoHeart**  
Automatic defibrillator

Art. AE801Z12  
**Wall cabinet**  
For TecnoHeart with front  
glass protection

## BUNDLE 8

**TecnoHeart +  
Wall cabinet**

List Price ~~1.670€~~

Offer Price

**1.400€**



Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp

Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit

Art. SB1A305+1ZB3A0006  
**Kit 50 disposable bags  
+ 5 Flovac connectors**

Art. 1498E  
**Kyri bridge**  
Variable height stainless  
steel bridge table


**BUNDLE 9**

**Orion 40DS +  
Kyri DSS +  
Kyri bridge =**

List Price ~~9.837€~~

Offer Price

**6.500€**

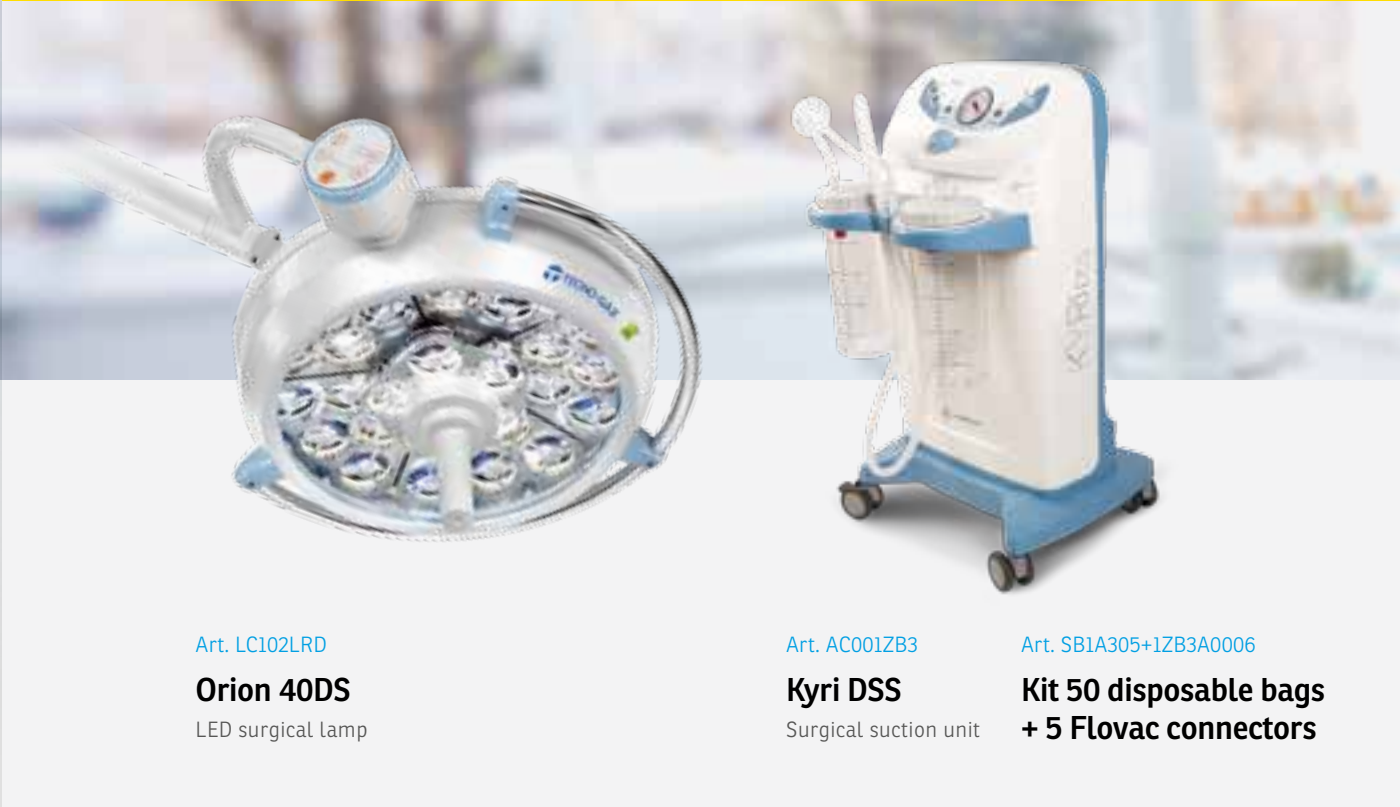

Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp

Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit

Art. SB1A305+1ZB3A0006  
**Kit 50 disposable bags  
+ 5 Flovac connectors**

**BUNDLE 10**

**Orion 40DS +  
Kyri DSS =**

List Price ~~9.157€~~

Offer Price

**6.000€**



Art. MO001STY  
**Smart Trolley**  
Electromedical

Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp

Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit

Art. SB1A305+1ZB3A0006  
**Kit 50 disposable bags  
+ 5 Flovac connectors**

Art. 1498E  
**Kyri bridge**  
Variable height stainless  
steel bridge table

**BUNDLE 11**

**Smart Trolley +  
Orion 40DS + Kyri  
DSS + Kyri bridge =**

List Price ~~11.355€~~

Offer Price

**7.900€**





**BUNDLE 12**

**Kyri + Orion40DS  
+ Kyri DSS + Kyri  
bridge + ErgoSeat**

List Price ~~21.111€~~

Offer Price

**14.500€**



Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp



Art. PC1..POC  
**Kyri**  
Unit for lengthy surgical  
and dental treatments



Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit

Art. SB1A305+1ZB3A0006  
**Kit 50 disposable  
bags + 5 Flovac  
connectors**



Art. 1498E  
**Kyri bridge**



Art. SD1..SED + SD999SED  
**ErgoSeat**  
Ergonomic work chair



Art. PC1..POC  
**Kyri**  
Unit for lengthy surgical  
and dental treatments



Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp



Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit

Art. SB1A305+1ZB3A0006  
**Kit 50 disposable bags  
+ 5 Flovac connectors**



Art. 1498E  
**Kyri bridge**



Art. M0006STY  
**Structure Parking x 4**



Art. SD1..SED + SD999SED  
**ErgoSeat**  
Ergonomic work chair

Art. M0001STY+M0002STY+M0014STY  
**Smart Trolley**  
3 Smart Trolleys,  
Electromedical, Intra Operatory and Postoperative



*Electromedical*



*Intraoperative*



*Postoperative*

**BUNDLE 13**

**Kyri + Orion 40DS +  
Kyri DSS + Kyri bridge +  
ErgoSeat + 3 Smart  
Trolleys + Structure  
Parking x 4 =**

List Price **~~27.173€~~**

Offer Price

**18.900€**





Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp

## BUNDLE 14

**Kyri + Kyri Cart +  
Orion 40DS + Kyri  
DSS + Kyri bridge +  
ErgoSeat =**

List Price

~~30.329€~~

Offer Price

**22.500€**



Art. PC1..POC  
**Kyri**  
Unit for lengthy surgical  
and dental treatments

Art. CA001K01+CA801K01+CA803K01  
+CA804K01+CA805K01  
**Kyri Cart**  
with syringe 3F, Module for Turbine,  
Micromotor Module NSK LED, Additional  
Scaler Module SATELEC w/handpiece



Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit

Art. SB1A305+1ZB3A0006  
**Kit 50 disposable bags  
+ 5 Flovac connectors**



Art. 1498E  
**Kyri bridge**



Art. SD1..SED + SD999SED  
**ErgoSeat**  
Ergonomic work chair



**BUNDLE 15**

**Kyri + Kyri Cart +  
Kyri DSS + Kyri  
bridge + Orion 40DS +  
ErgoSeat + 2 Smart  
Trolleys + Structure  
Parking x4 =**

List Price **~~34.873€~~**

Offer Price

**24.900€**

Art. PC1..POC

**Kyri**

Unit for lengthy surgical  
and dental treatments

Art. LC102LRD

**Orion 40DS**

LED surgical lamp

Art. CA001K01+CA801K01+CA803K01  
+CA804K01+CA805K01

**Kyri Cart**

with syringe 3F, Module for Turbine,  
Micromotor Module NSK LED, Additional  
Scaler Module SATELEC w/handpiece

Art. AC001ZB3

**Kyri DSS**

Surgical suction unit

Art. SB1A305+IZB3A0006

**Kit 50 disposable bags  
+ 5 Flovac connectors**

Art. M0006STY

**Structure Parking x 4**

Art. M0002STY+M0014STY

**Smart Trolley**

2 Smart Trolleys,  
Intra Operatory and Postoperative

*Intraoperative*

*Postoperative*

Art. SD1..SED + SD999SED

**ErgoSeat**

Ergonomic work chair

Art. 1498E

**Kyri bridge**





Art. LC102LRD  
**Orion 40DS**  
LED surgical lamp



Art. PCL..POC  
**Kyri**  
Unit for lengthy surgical  
and dental treatments



Art. SD1..SED + SD999SED  
**ErgoSeat**  
Ergonomic work chair



Art. 1498E  
**Kyri bridge**



Art. AC001ZB3  
**Kyri DSS**  
Surgical suction unit



Art. CA001K01+CA801K01+  
CA803K01+CA804K01+CA805K01  
**Kyri Cart**  
with syringe 3F, Module for  
Turbine, Micromotor Module NSK  
LED, Additional Scaler Module  
SATELEC w/handpiece



Art. GI001K02  
**Kyri Water  
Unit**

## BUNDLE 16

**Kyri + Orion40DS + Kyri  
DSS + Kyri Cart + Kyri  
Water Unit + Kyri bridge +  
ErgoSeat**

List Price

~~33.716€~~

Offer Price

**24.900€**



*Anticipate, propose and guide the change.*

# WE MANUFACTURE THE FUTURE

*BUILDING THE FUTURE,  
DAY AFTER DAY,  
IDEA AFTER IDEA.*



HEADQUARTERS



FACTORY 2



AUTOCLAVE PRODUCTION LINE



**Tecno- Gaz S.p.A.**

Strada Cavalli, 4 - 43038 - Sala Baganza - Parma - Italy

Ph. +39 0521 83.80 Fax +39 0521 83.33.91 - [www.tecnogaz.com](http://www.tecnogaz.com)

Cap. Soc. € 280.000 i.v. C.F. e P.IVA/VAT IT00570950345 - R.E.A. PR 138927 Iscr. Reg. Impr. PR 10061

All rights reserved by Tecno- Gaz S.p.A. Changes in images or content may be made without prior notice. Tecno- Gaz S.p.A. is not liable for damage resulting from incorrect or incomplete information reported in this document. The images used in this document are for illustrative purposes only.

Visit our Website  
[www.tecnogaz.com](http://www.tecnogaz.com)

For More Information:  
**00 39 0521.8380**

